



LABORATORY LINES

AUTOCLAVE



Kalstein®

A different accompaniment, at your service





Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ❤️

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



Digital Display Vertical Type Pressure Steam Autoclave

Model Series YR05676

Application

YR series vertical pressure steam autoclave is the equipment that disinfects and sterilizes the objects with saturated steam, its structure is horizontal and heating way is electric heating. It is composed of sterilization room, shell, door, steam generator, pipe system, control system, digital display system and power-supply system.

It is widely used in CSSD, lab, scientific research institution, pharmaceutical, food industry and is suitable for items which are high temperature-resisting and high pressure-resisting including unwrapped or wrapped instrument, dressing, rubber, liquid, plant-in instrument, porous item, hollow instrument to be sterilized.

Product Features

- ✓ Fully stainless steel 304 structure.
- ✓ Wing nut locking mechanism.
- ✓ Pressure gauge, real-time display of pressure.
- ✓ Microprocessor control, digital display of working status: pressure, time, temperature, error code, touch type key.
- ✓ Auto discharge the cool air, and steam discharging automatically after sterilizing.
- ✓ Self-inflating type seal.
- ✓ Automatically shut off with beep reminding after sterilization.
- ✓ Standard PT/TT test interface is provided for convenient testing.
- ✓ With two stainless steel sterilization baskets.
- ✓ With draining output tap and moving wheels.

Safety Protection Device

- ✓ Equipped with water level sensor, alarm and cut-off if water lacking.
- ✓ The program will be stopped if water level lower than the low limit.
- ✓ The buzzer alarms when a cycle of operation is completed.
- ✓ Equipped with over-temperature and over-pressure auto-protection.
- ✓ Over current protector.
- ✓ Door safety interlock device.
- ✓ Safety valve for automatically released.
- ✓ Emergency exhaust valve.



Model	Series YR05676
Control	Microprocessor
Display	Digital display
Material	Fully stainless steel 304
Working pressure	0.22Mpa
Working temperature	134°C
Working temperature range	105 - 134°C
Sterilization time range	0 - 99 min or 0-99 hour 59min
Temperature accuracy	0.1°C
Heat average	≤1°C
Voltage	110-240V. 50/60Hz
Steam generator	Stainless steel 304
Loading	2 pcs basket
Caster	Casters with brake



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S.



Digital Display Vertical Type Pressure Steam Autoclave

Model Series YR05676



Model	Chamber size (Dia×D)(mm)	Outer size (mm)	Packing size (mm)	G.W. / N.W.	Power
YR05676	35L / Ø318×450 mm	450×450×1010 mm	570×570×1150 mm	72 kg / 56 kg	2.5 KW
YR05677	50L / Ø340×550 mm	510×470×1130 mm	590×590×1280 mm	88 kg / 68 kg	3.0 KW
YR05678	75L / Ø400×600 mm	560×560×1130 mm	650×650×1280 mm	100 kg / 80 kg	4.5 KW
YR05679	100L / Ø440×650 mm	540×560×1250 mm	680×680×1370 mm	110 kg / 85 kg	4.5 KW
YR05680	120L / Ø480×660 mm	600×640×1140 mm	730×730×1370 mm	125kg / 100kg	6.0 KW
YR05681	150L / Ø510×740 mm	670×690×1130 mm	760×760×1370 mm	140kg / 110kg	6.0 KW



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

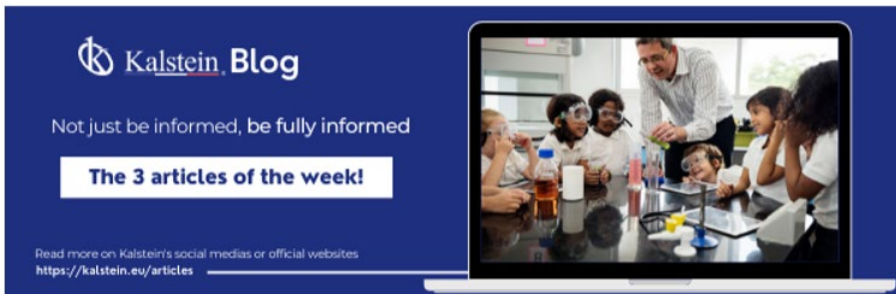

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

refrigerator is a device that is one of the most common pieces of equipment used in laboratories. It is a controlled environment (refrigerated space), various fluids and substances

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?


Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago

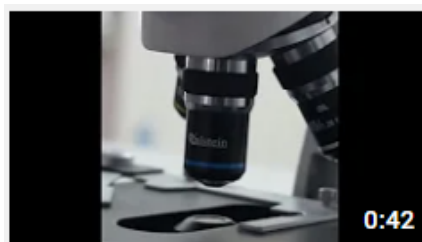

Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

