




LABORATORY SECTOR

PH METERS



Kalstein®

A different accompaniment, at your service

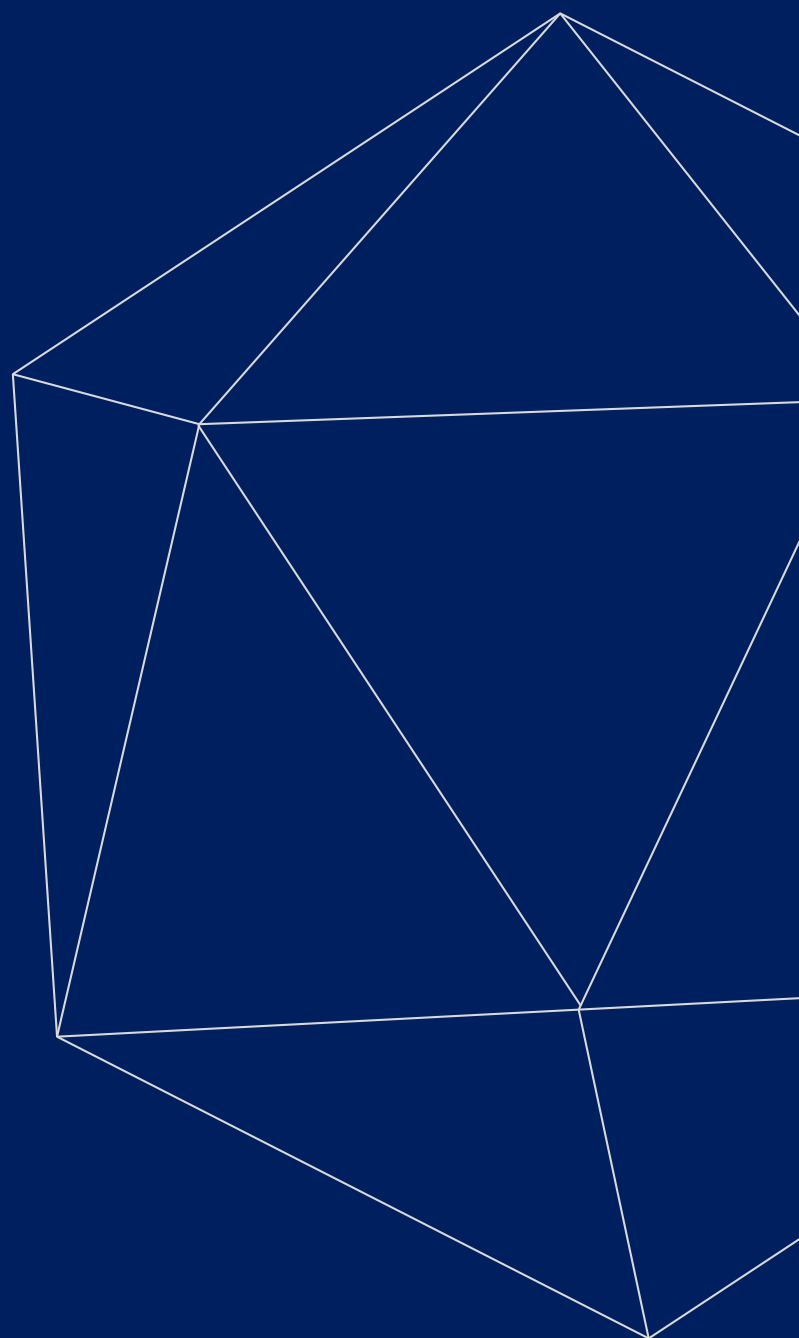



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® PH Meters
Benchtop Dissolved Oxygen Meter
Model YR01832



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® PH Meters
Benchtop Dissolved Oxygen Meter
Model YR01832




Product Description:

Professional benchtop dissolved oxygen meter, 7 inch TFT display, stores up to 1000 data sets, setup menu contains 17 options, accuracy: 0.2mg/L.



Features:

- ✓ Professional benchtop dissolved oxygen meter is equipped with a 7 inch TFT display.
- ✓ 1 or 2 points calibration using the air-saturated water or zero oxygen solution.
- ✓ Salinity and barometric pressure compensations eliminate the measurement error.
- ✓ Auto-Read function senses and locks the measurement endpoint.
- ✓ Timed Interval Readings sends the measurement data to a PC or printer.
- ✓ Limit Alarm automatically alerts when reading exceeds the specified range.
- ✓ Calibration due alarm prompts user to calibrate the meter regularly.
- ✓ Calibration report shows the date, time, calibration point, etc.
- ✓ Password protection prevents the unauthorized calibration and settings.
- ✓ Expanded memory stores or recalls up to 1000 data sets.
- ✓ Reset function automatically resumes all settings back to the factory defaults.

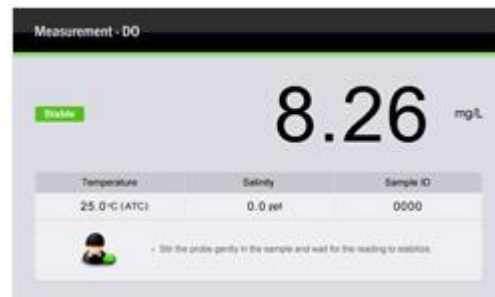

#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® PH Meters

Benchtop Dissolved Oxygen Meter

Model YR01832



Technical Description:

Model	YR01832
Dissolved Oxygen	
Range	0.00 to 20.00mg/L
Resolution	0.01mg/L
Accuracy	±0.2mg/L
% Saturation of Oxygen	
Range	0.0 to 200.0%
Resolution	0.1%
Accuracy	±2.0%
General Specifications	
Calibration Points	1 or 2 points
Temperature Compensation	0 to 50°C/32 to 122°F, automatic
Pressure Correction	60.0 to 112.5kPa/450 to 850mmHg, manual
Salinity Correction	0.0 to 50.0g/L, manual
Measurement Modes	DO
Stability Criteria	Standard or high-accuracy
Calibration Due Alarm	1 to 31 days or off
Interval Readings	10, 30, 60 seconds, 10, 30 minutes or off
Password Protection	4 digits
Memory	Stores up to 1000 data sets
Output	USB communication interface
Connector	6-pin nemi-DIN
Display	7" TFT LCD
Power	DC12V, using AC adapter, 220V/50Hz
Dimensions	240(L)×220(W)×80(H)mm
Weight	1.7kg
Order Code	Accessories
YR01832	DO100 dissolved oxygen electrode
	DO electrolyte solution
	DO membrane cap
	DC12V power adapter



#Letsgivemore 

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

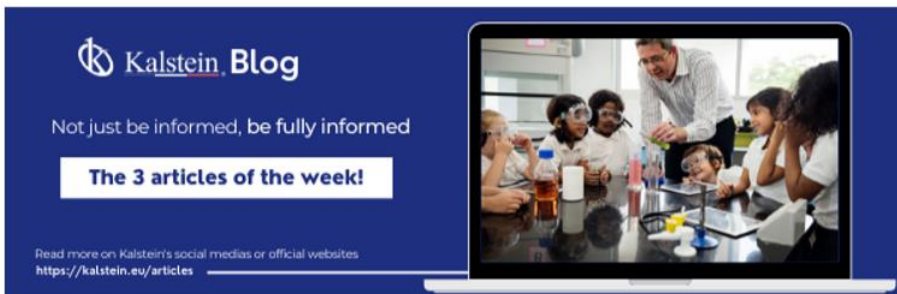

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89

Cel: +33 (0) 6 80 76 07 10

sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

