




MEDICAL LINES

ELECTRIC MEDICAL BED



Kalstein®

A different accompaniment, at your service

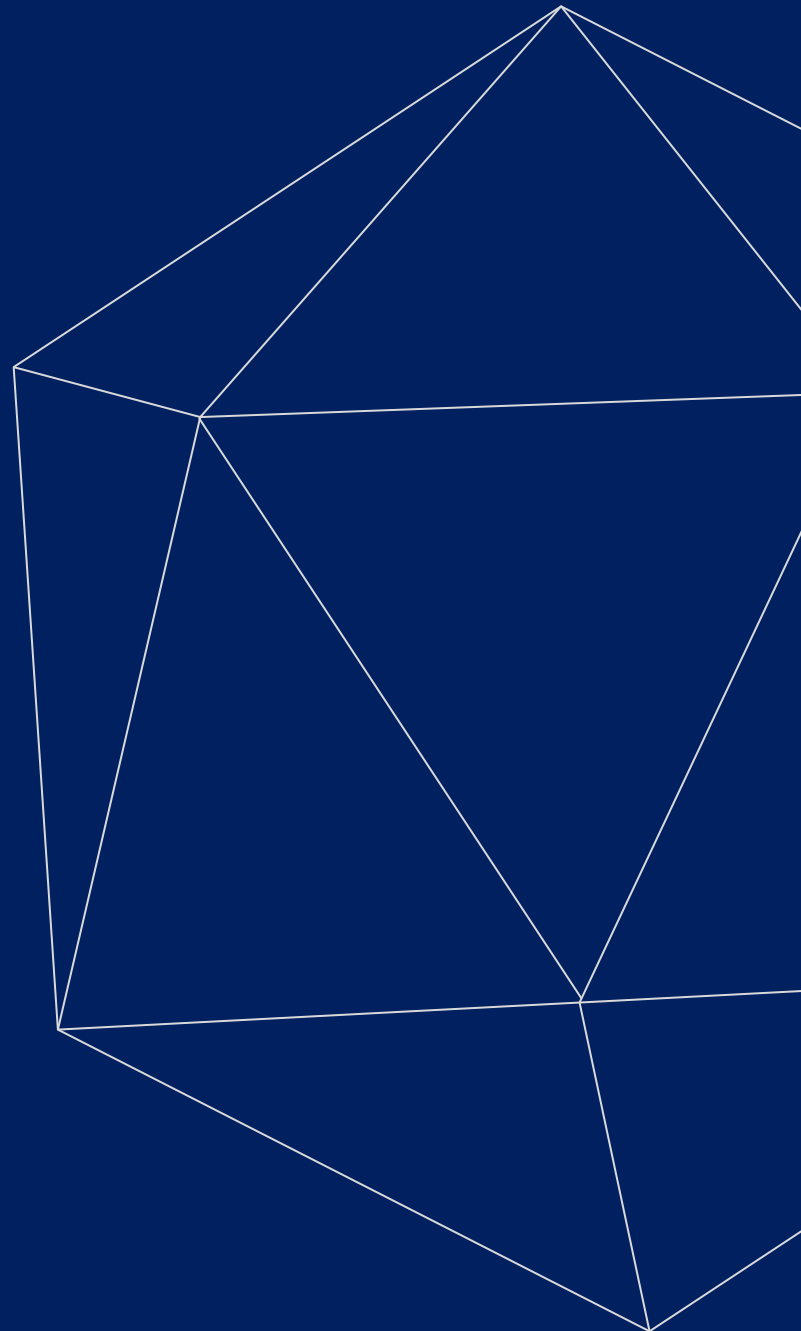



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



Medical Electric Hospital ICU Bed Model YR68c




ELECTRIC SYSTEM

American motor, with UL/ROHS/EN approval.
Noise under 50dB, with lifetime 20,000times movement.

UL / ROHS / EN

50 | 20,000_{times}

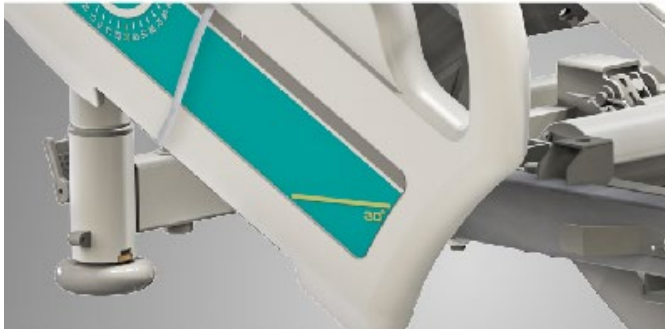

BACK REST 30DEGREE POSITION

Help nurse to observe the bed position from distance,
many research shows 30degree position will max.
reduce the chance of lung infection.

FUNCTION PROFESSION



Back-rest




Knee-rest



Height Adjustment



Bed extension

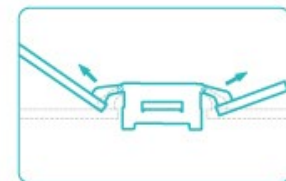

Auto-regression



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S




Medical Electric Hospital ICU Bed Model YR68c



Important Quality Control

- * Follow IEC60601 Safety Standard
- * American motor, with UL/ROHS/EN approval. Noise under 50dB, with lifetime 20,000times movement.
- * Motor protection:
 - * The elevators are equipped with electric protection against overload, therefore the bed will automatically stop when an elevator is overloaded.
 - * Once the overload is removed, the bed can continue to work without interference.
 - * One-button emergency stop switch can quickly stop all going movement
 - * Use the GO button to activate the bed; protects the patient against an injury caused by any unintentional movement of the bed. Control system is automatically locked after 3 min
 - * The designed guardrail fully considered to provide enough space and hand grasp for patient to get out of bed, and provide extremely high safety.
 - * Two steps unlocking system, to prevent patients from accidentally lowering the guardrail (patented design).
- * NURSE CONTROLS WITHIN REACH
- * Clear and easy to use controllers which is at the patient's disposal
- * Auto-Regression 10+6 cm
- * When repositioning the backrest and thigh rest, the Y6Y8C makes the back plate retreat 10cm and the thigh plate retreats 6cm so to minimises the pressure and shear forces in the patient's pelvic region. This results in a significant reduction in the development of pressure ulcers.
- * Respects easy patient mobilisation:
 - * One-button for patient mobilisation position -back plate raised, double hole design on back part guardrail, which supports patient self standing
 - * CPR emergency reset handle to immediately flatten the bed back part; located on both sides of the bed.
 - * Electric CPR button is settled inside control panel and nurse control panel, one button to make the back plate flat.

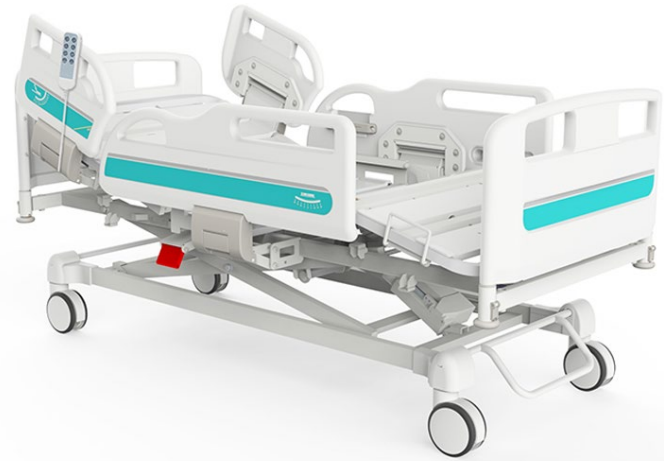

- * One button nursing operating position: 800mm, helping the nurse to give fast treatment in a very comfortable height
- * The lowest position indicator design allows the medical staff to observe the bed safe position during night or long distance; green light indicate lowest bed height
- * Backrest adjustment
- * Thigh rest adjustment
- * Auto-contour
- * Height adjustment
- * Bed extension
- * X-Ray translucent bed surface
- * Auto-regression
- * Cardiopulmonary resuscitation
- * Controls in the siderails: Numerous studies have shown that the built-in controller located in the head guardrail is easier to reach, allowing the patient to adjust the basic position of the bed through the controller. (optional)

Parameter

Bed dimentions	L2130*W1050*H350-800mm
Bed platform	L1925*W900mm
Back rest	75°±10°
Leg rest	35°±10°
Packing size	2120*1080*450mm
Mattress size	L1925*880*10mm
Bumper wheel	4xΦ80mm



#Letsgivemore ♥


Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tit: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Medical Electric Hospital ICU Bed Model YR68c

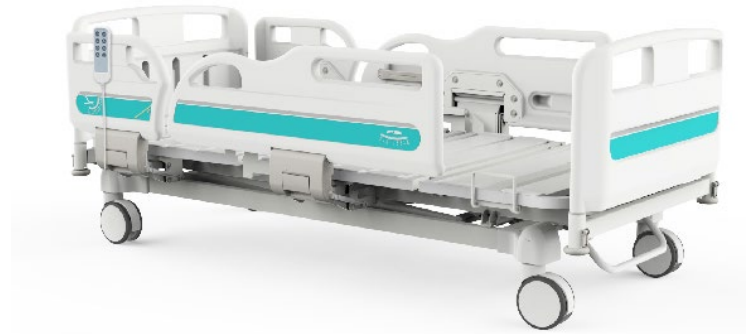

- * Handheld controller for patients and nurses, ergonomic design. Built-in LED torch, backlit keyboard design for night use
- * Satellite control panel for nurse: flexible spiral cable for easy mobility and operation in optimal position. The nurse can control all standard functions and automatic chair position. This panel is an optional accessory.
- * Patented design Detachable head & foot board, using the eccentric locking structure, one step quick removable
- * Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall
- * Patent angle indicator, clearly show back-rest
- * extended bed platform up to 20CM, comfortable for the tallest of patients (option)
- * JOINT-FREE TECHNOLOGY
- * Removable plastic covers formed from 9 removable blow moulded lamella to minimise the slots where unwanted infections may hide.
- * 5th and 6th extra side rail, for 100% fully protection(option)
- * Castors with TPR tire, no worn-out after running 30KM, united forming without bolts
- * Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time
- * Plastic drainage hook, movable
- * Full SS IV pole, bearing 5kg
- * Socket on side rail for placing the patient cell phone or glasses
- * Back rest 30degree position - help nurse to observe the bed position from distance, many research shows 30degree position will max. reduce the chance of lung infection.
- * Mattress
- * Ergoframe® against pressure ulcers
- * Easy mobilisation
- * Available with a number of controls
- * Soft drop function
- * Easy clean concept
- * Safety features such as the Extender or Protector
- * Low TCO – total cost of ownership
- * Universal bed for a low price
- * Overload protection
- * Baosteel® from Fortune Global 500 company
- * Environmental plastic can be sterilized in 100°C, tensile up to 30MPa
- * 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact
- * Panasonic Robotic ensure
- * 360° full smooth welding



#Letsgivingmore ❤️

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



EASY PATIENT MOBILISATION

Respects easy patient mobilisation:
One-button for patient mobilisation position
-back plate raised, double hole design on
back part guardrail, which supports patient
self standing.




SAFETY GUARDRAIL

The designed guardrail fully considered to provide
enough space and hand grasp for patient to get out
of bed, and provide extremely high safety.
Two steps unlocking system, to prevent patients from
accidentally lowering the guardrail (patented
design).

AUTO-REGRESSION

Auto-Regression 10+6 cm
When repositioning the backrest and thigh rest,
makes the back plate retreat 10cm and
the thigh plate retreats 6cm so to minimises the
pressure and shear forces in the patient's pelvic
region. This results in a significant reduction in
the development of pressure ulcers.

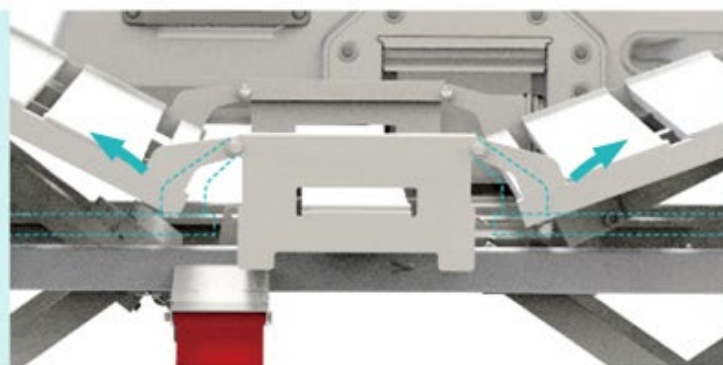


SAFETY

Follow IEC60601 Safty Standard
One-button emergency stop switch can
quickly stop all going movement.

CPR

CPR emergency reset handle to immediately
flatten the bed back part; located on both sides
of the bed.
Electric CPR button is setted inside control
panel, one button to make the back plate flat.

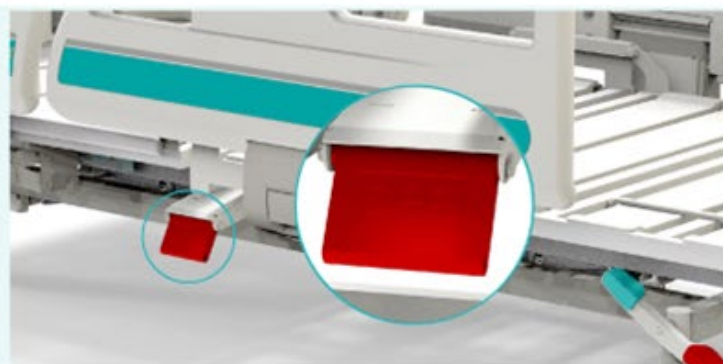

#Letsgivemore

Thanks to your purchase, a donation will be made
to a non-profit foundation that fights against breast
cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



RICH PARTS



Castors with TPR tire, no worn-out after running 30KM, united forming without bolts

Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 times.



JOINT-FREE TECHNOLOGY

Removable plastic covers formed from 9 removable blow moulded lamella to minimise the slots where unwanted infections may hide.



Patent angle indicator, clearly show back-rest and trolley position



Patented design Detachable head & foot board, using the eccentric locking structure, one step quick removable



Extended bed platform up to 20CM, comfortable for the tallest of patients (option)



Mattress Holder.
Waterproof, antibacterial mattress.



Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall



5th and 6th extra side rail, for 100% full protection



Socket on side rail for placing the patient's cellphone or glasses



Plastic drainage hook, movable



Full SS IV pole, bearing 5kg



#Let'sgivemore 


Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710 / +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S



Articles Kalstein
(Social Networks and Official Websites)

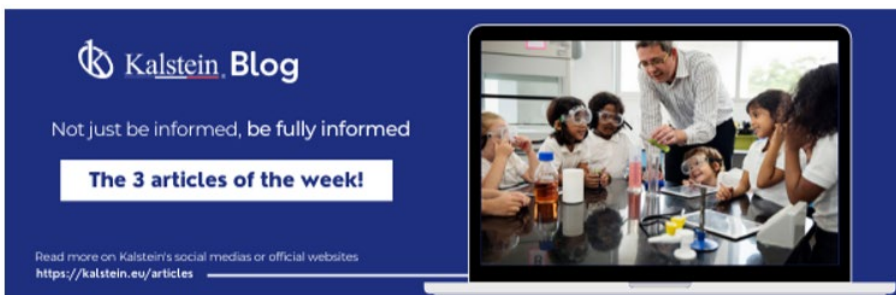

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
A microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form...

Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
A refrigerator is a device that is one of the most common pieces of equipment used in laboratories to maintain, in a controlled environment (refrigerated space), various fluids and substances...

Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



MONTPELLIER - FRANCE
HEADQUARTER
Optimum Business Center,
450 Rue Baden Powell,
34000 Montpellier - France.

Fax: +33 (0) 4 67 15 88 49
Tif: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

Official Websites for Latin America.

