




LABORATORY SECTOR

HOTPLATES AND HOTPLATE STIRRERS



Kalstein®

A different accompaniment, at your service

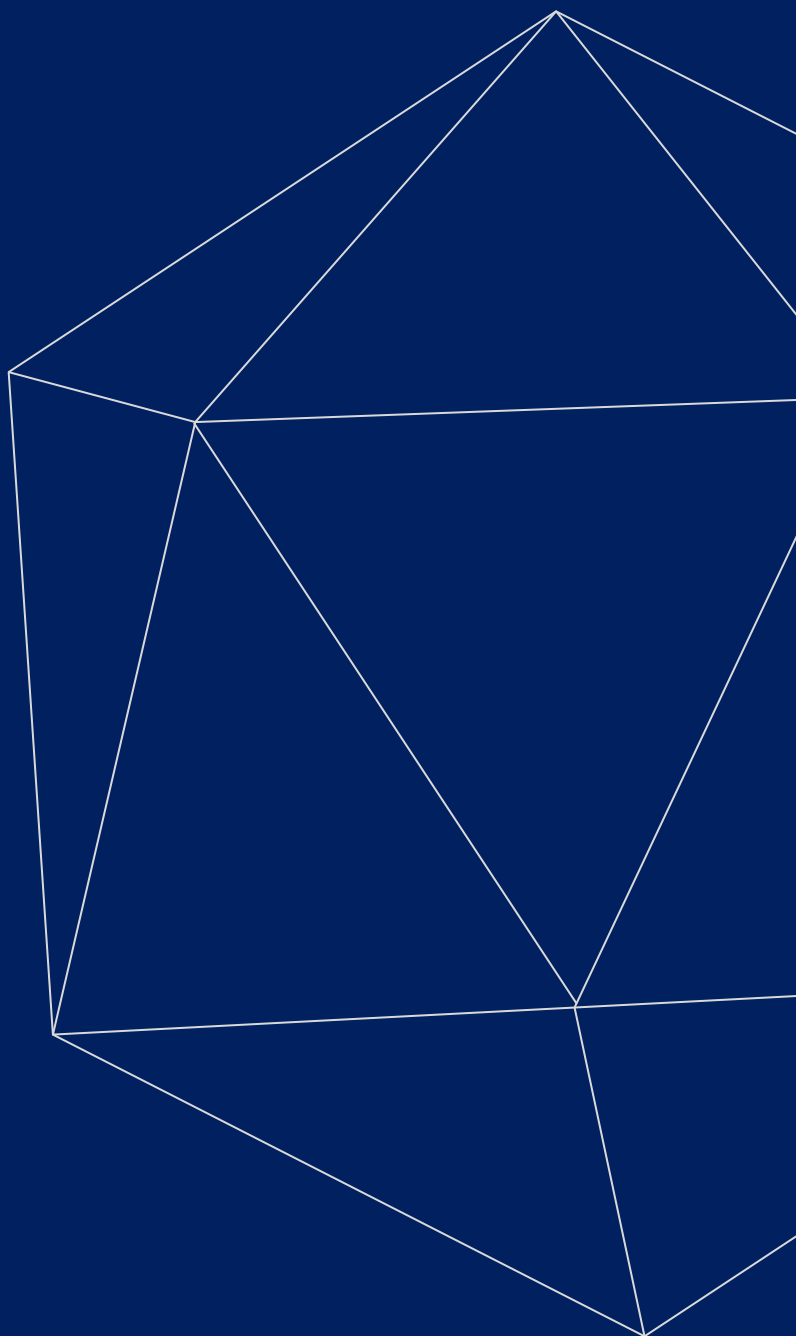



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Hotplates and Hotplate Stirrers
Hotplate Magnetic Stirrer Bar With Heating
Model YR02943



#Letsgivingmore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Hotplates and Hotplate Stirrers

Hotplate Magnetic Stirrer Bar With Heating

Model YR02943



Product Description:

Magnetic stirrers are often used in chemistry and biology. They are preferred over gear-driven motorized stirrers because they are quieter, more efficient, and have no moving external parts to break or wear out (other than the simple bar magnet itself). Because of its small size, a stirring bar is more easily cleaned and sterilized than other stirring devices. They do not require lubricants which could contaminate the reaction vessel and the product. They can be used inside hermetically closed vessels or systems, without the need for complicated rotary seals. Magnetic stirrers may also include a hot plate or some other means for heating the liquid.

Product Features:

- ✓ Very compact and simple to operate
- ✓ Digital Display temperature and speed, Heating temperature: Rt-380 Degree and mixing speed :100-1800 RPM
- ✓ Heat and corrosion resistant aluminum alloy housing for durability
- ✓ Ceramic coated aluminum heating plate for great heat distribution and corrosion resistance
- ✓ Specially designed housing prevents spill from damaging the electronic components inside
- ✓ Secondary heat insulation plate prevents heat loss and eliminates heat damage to electronic components
- ✓ Adjustable temperature and rotation speed



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Hotplates and Hotplate Stirrers

Hotplate Magnetic Stirrer Bar With Heating

Model YR02943



Technical Description:

Model	YR02943
Volume	10L
Heating Plate (MM)	230*230
Heating Power(W)	750
Heating Plate Material	Ceramic Coated Aluminum (White or Black)
Temperature Range	Rt-380°C
Mixing Speed	100-1800RPM
Temperature Control	PID Auto Control
Temperature Display (°C)	Digital LCD Display
Temperature Accuracy	+/- 1°C
Thermocouple Probe(Sensor)	PT100, with Inner and Outer Thermocouple
Permissible Working Temperature	Rt to 40°C
Voltage	220V/ 110V Optional
Dimensions (mm)	500 x 400 x 350
Weight (Kg)	7.5



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:

 **Fundación
Jacinto Convit**
<https://www.jacintoconvit.org/>




<https://onetreepanted.org>

 **Humatem**
Solidarité et coopération
biomédicale
<https://www.humatem.org/>


**MANIAPURE
FOUNDATION USA**
<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

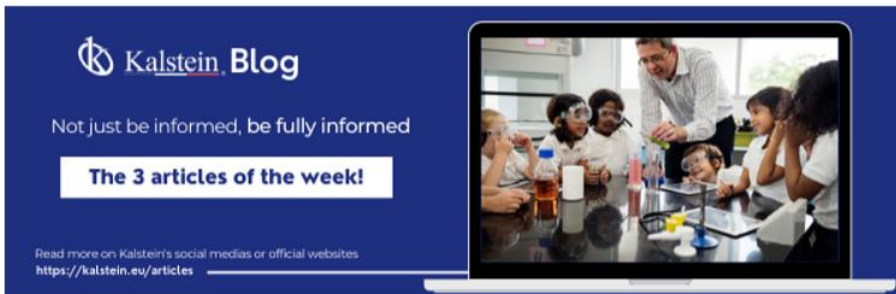

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
microscopes in a laboratory are one of the most sensitive pieces of equipment that require to maintain in a controlled environment (refrigerated space), without fluctuations in humidity and temperature.
Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
Refrigerator in a laboratory is one of the most sensitive pieces of equipment that require to maintain in a controlled environment (refrigerated space), without fluctuations in humidity and temperature.
Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER

2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

