




MEDICAL LINES

INFANT PHOTOTHERAPY INCUBATOR



Kalstein®

A different accompaniment, at your service

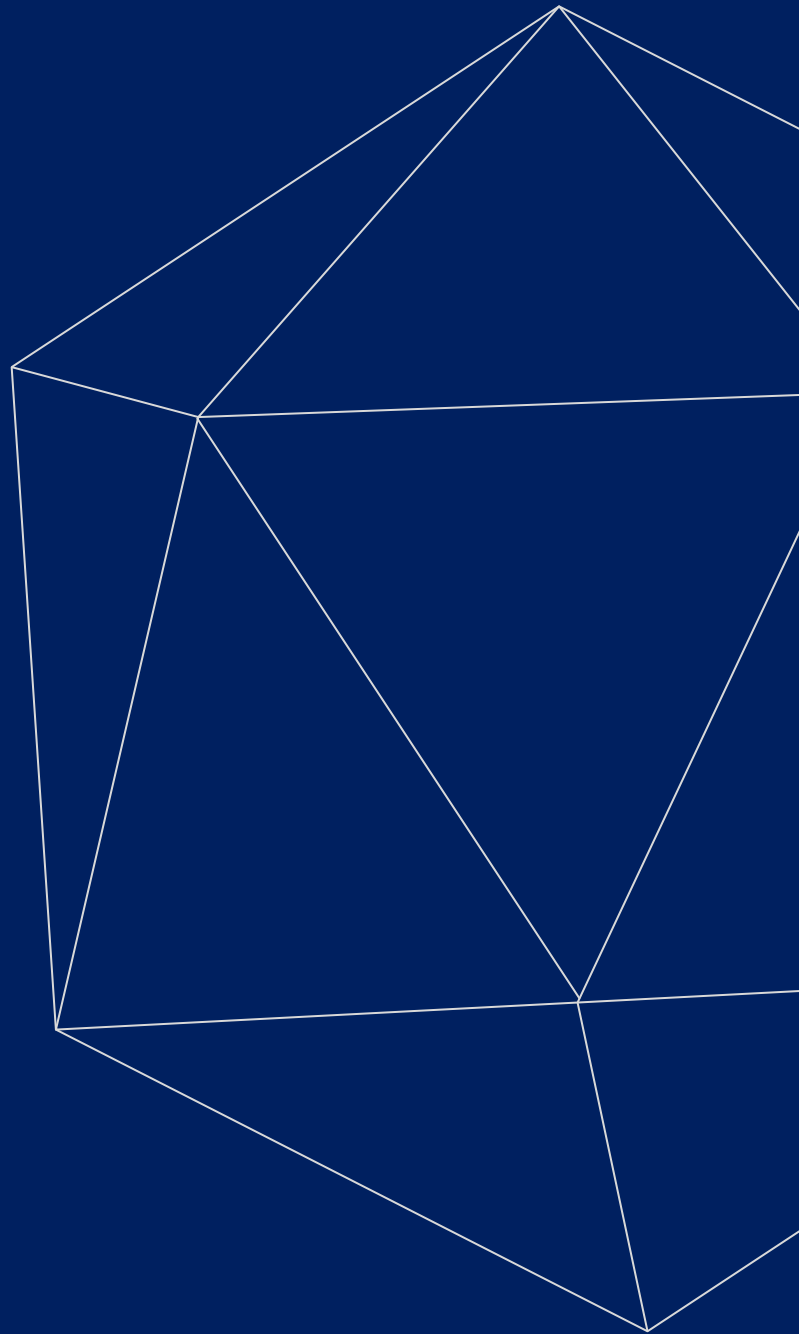



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Infant Phototherapy Incubator

YR02183



Kalstein's our medical use Baby incubator for jaundice Model YR02183 keeps the neonate patients in stable environment and offers best jaundice therapy in neonatal intensive care unit. The design of phototherapy incubator gives the fragile patients the best comfortable environment. Guarantee an easy and friendly way for doctors and nurses.

Application Industry

We designed the infant phototherapy incubator is for fragile patients, and it has your exceptional requirements and best care solutions. Widely applicable in all hospitals and clinics and optimize for preterm births or for some ill full-term babies or jaundiced baby.

Features

- Phototherapy light intensity is constantly adjustable
- Time of bilirubin phototherapy can be paused, recorded and accumulated
- Phototherapy treatment duration can be pre-set and it automatically shuts off when a pre-set duration is reached
- Working life of LED phototherapy unit is lasting 20,000 hours
- Silicone mattress protect neonate from cauma while phototherapy treatment
- Microprocessor servo-controller ensures a precise and reliable temperature control
- Set temperature, monitored air temperature and baby skin temperature and heating power are distinguished in display
- Servo-controlled humidification system maintains the relative humidity at a constant level in the incubator
- Set and monitored humidity concentration are displayed separately.
- Cupboard cabinet.



Features

- Ultrasonic and low temperature nebulized humidification protects neonate from bacterias pollution
- > 37°C override temperature function
- Self-checking function for safety affirmative
- Audible and visible alarm indications for multiple failures with 3 priorities
- Chord alarm notes nurse farther than 8meters, but leaves asleep baby uninterrupted
- Aluminous framework for constant temperature keeping
- Continuous adjustment of mattress inclination up to 10°
- Drawer-type water reservoir is easily removable for cleaning and disinfection
- Nature air circulation system



Air and baby skin control
Humidity control



#Let'sgivemore ❤️

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S.



KALSTEIN® Infant Phototherapy Incubator

YR02183



Ultrasonic and low temperature nebulized humidification for YR02182 - YR02183 - YR02184



Vertical height adjustment with ergonomic design for YR02184 only



X type base for YR02182 - YR02183



Basic base for YR02182 only



High illumination intensity treatment with double-sided phototherapy for YR02182 - YR02183 - YR02184




User friendly control interface



Baby bed can be adjusted up to 12° for YR02182 - YR02183 - YR02184



Radiolucent silicone mattress (optional) for YR02182 - YR02183 - YR02184

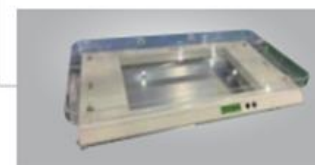

Baby scale bed (optional) for YR02182 - YR02183 - YR02184

YR02183



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
2 Rue Jean Lantier • 75001 Paris •
+33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
KALSTEIN FRANCE, S. A. S



KALSTEIN® Infant Phototherapy Incubator

YR02183



Technical Specifications

Model	YR02183
Air temperature control	●
Control of skin temperature	●
37 ° air temperature control	●
The baby bed can be tilted	●
Humidity display	●
Humidity control	●
Superior Phototherapy Unit	●
Lower Phototherapy Unit	●
Side door	○
Aluminum tank for heat preservation	●
Windows and access ports operation	1 port, 4 windows
Cabinet	●
Adjustable working height	○
Additional rear door	○
Baby scale bed	○
Silent DC motor	○

Optional Accessories

- Side door, front-rear door
- Double wall hood
- Vertical height adjustable
- Infant phototherapy unit (OSRAM tubes)



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
 2 Rue Jean Lantier • 75001 Paris •
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •
<https://kalstein.eu>
 KALSTEIN FRANCE, S. A. S





A different accompaniment, at your service

WITH THE ACQUISITION OF A KALSTEIN EQUIPMENT

YOU MAKE A CONTRIBUTION TO:

#LETSGIVEMORE



Fundación Jacinto Convit

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération biomédicale

<https://www.humatem.org/>




MANIAPURE FOUNDATION USA

<https://www.maniapurefoundation.org/>

Thank you!

Articles Kalstein
(Social Networks and Official Websites)

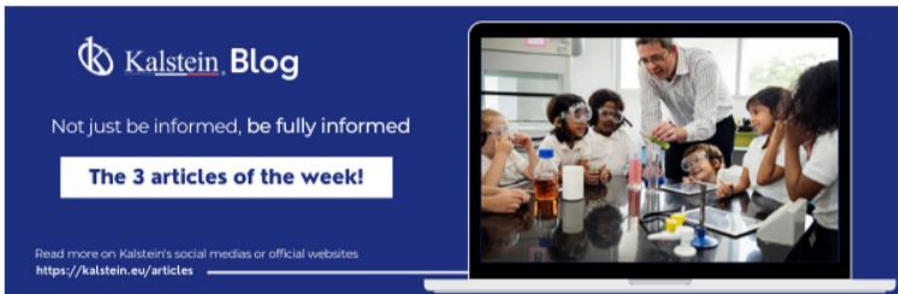

Kalstein Blog
Not just be informed, be fully informed
The 3 articles of the week!
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



1
Read this ARTICLE NOW
How does the climate affect your microscope?
microscopes in a laboratory are one of the most sensitive pieces of equipment that require to maintain in a controlled environment (refrigerated space), without fluctuations in temperature and humidity.
Kalstein Blog
Article of the week:
How does the climate affect your microscope?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW
What is the importance of the laboratory refrigerator?
Refrigerator in a laboratory is one of the most sensitive pieces of equipment that require to maintain in a controlled environment (refrigerated space), without fluctuations in temperature and humidity.
Kalstein Blog
Article of the week:
What is the importance of the laboratory refrigerator?
Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances
23 views • 3 years ago



Kalstein's Microscopes
445 views • 1 week ago



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>

