




# LABORATORY SECTOR

REFRIGERATORS AND FREEZER



**Kalstein**®

A different accompaniment, at your service

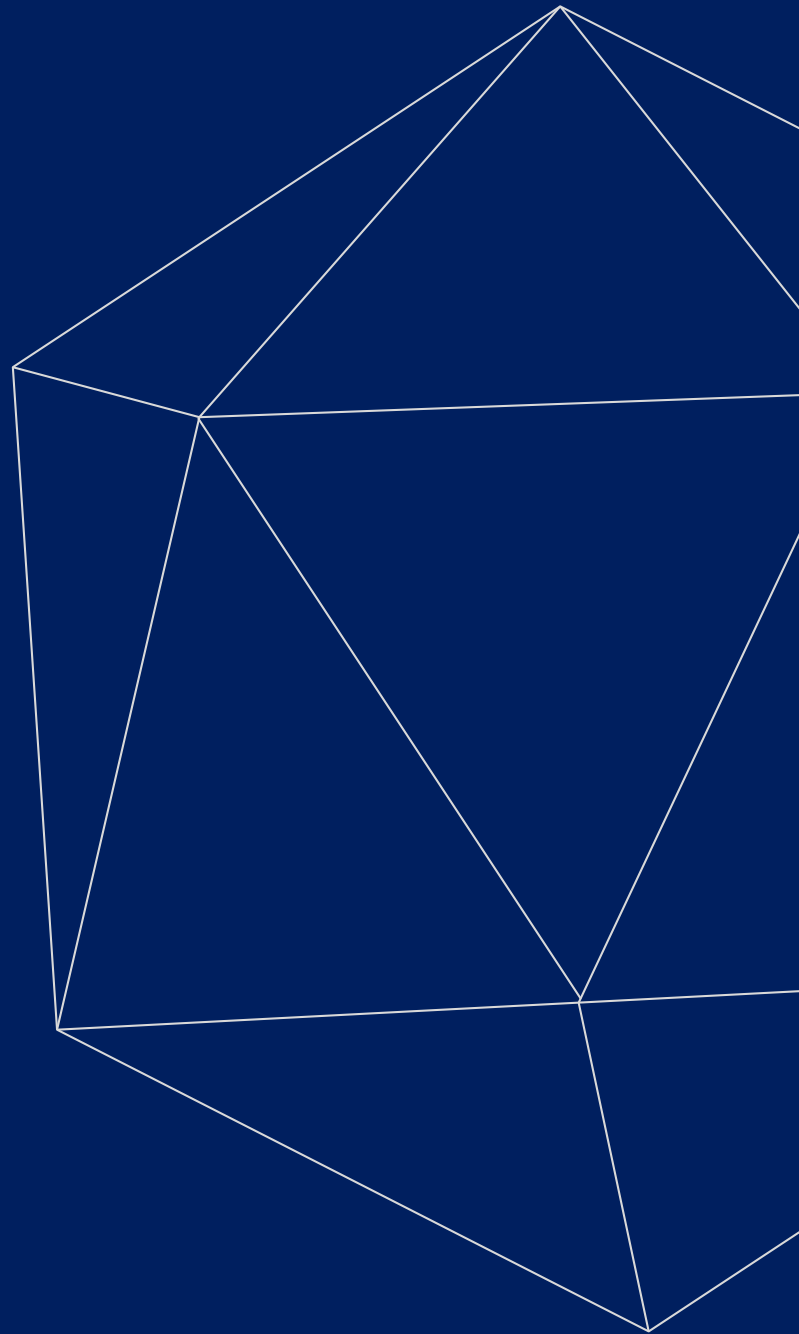



# Kalstein®

A different accompaniment, at your service

OUR SERVICES

## Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



### Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



### Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



### #Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



### Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



### Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



### Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



# KALSTEIN® Refrigerators and Freezer Medical Refrigerator Model YR05083



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
KALSTEIN FRANCE, S. A. S



# KALSTEIN® Refrigerators and Freezer

## Medical Refrigerator

### Model YR05083



#### Product Description:

The pharmacy refrigerators are suitable for installation in drug and vaccine stores, pharmaceutical companies, hospitals, epidemic prevention centers and clinics.

#### Hot Tips:

- ✓ Reliable control system
- ✓ Inside temperature 2~8°C
- ✓ Forced-air cooling system, Frost-free
- ✓ Service-free drainage system
- ✓ Wide voltage range, 187V~242V available

#### Refrigeration System:

- ✓ Imported brand compressor work with CFC-Free environmental refrigerant.
- ✓ High density insulation foam layer with better heat preservation effect.

#### Features:

- ✓ Adjustable temperature range 2-8°C
- ✓ Microprocessor control, adjustable temperature
- ✓ LED digital display and increment at 0.1 °C.
- ✓ Temperature variation within ±3°C
- ✓ Permanently lubricated cooling fan for safety and longevity
- ✓ High efficiency specialty compressor with known field reliability
- ✓ Forced- air cooling with optimized air distribution system designed
- ✓ To achieve maximum temperature uniformity and stability
- ✓ Optimized refrigeration system design for more effective cooling and speedy recovery

#### Humanized Design:

- ✓ Optimized space utilization design with multi – level and adjustable shelf height to accommodate storage requirement of pharmaceutical products of various kinds
- ✓ Standard with LED Interior light
- ✓ Safety lock to prevent unauthorized access requirement of pharmaceutical products of various kinds

#### Safety Guarantee:

- ✓ Build-in back-up battery to display temperature and to operate audible and visual alarm system for up to 8 hours without AC power
- ✓ Visual and audible alarm system: High/low temperature, thermostat failure alarm, power failure alarm, door ajar.
- ✓ Power failure protection: turn on delay of the cooling system when power failure, restart delay protection.
- ✓ Keyboard lock, password protection to avoid set parameter randomly.
- ✓ Wide available voltage range: 187V~242V.



#Letsgivemore ♥

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
KALSTEIN FRANCE, S. A. S



# KALSTEIN® Refrigerators and Freezer

## Medical Refrigerator

### Model YR05083



#### Technical Description:

Model	YR05083	
<b>Specificaitons</b>	Shelf/drawer	Shelves/2
	Cooling method	Forced Air Cooling System
	Defrost	Auto
	Refrigerant	R600a
	Power consumption(kWh/24h)	0.84
	Noise(db)	50
	Ambient temperature(°C)	10~32°C
	Temperature range(°C)	2~8°C
<b>Cooling system</b>	Compressor brand/QTY	China Brand/1
	Sensor Type	NTC
	Temperature controller	Microprocessor
	Display	Digital display
<b>Power</b>	Voltage/frequency(V/Hz)	220/50Hz
	Power(W)	72W
	Current(A)	0.59A
<b>Material</b>	Inside material	Color sprayed steel
	Outside material	Color sprayed steel
	Insulation	PURF
<b>Size</b>	Capacity(L/cu.ft)	60/2.11
	NT./GT.(kg)	33/37
	Exterior size(W*D*H)(mm)	495*450*710
	Interior size(W*D*H)(mm)	410*315*560
	Package size(W*D*H)(mm)	558*510*760
	CBM	0.22
		20GP/40GP/40HQ



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
2 Rue Jean Lantier • 75001 Paris •  
+33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
KALSTEIN FRANCE, S. A. S



**KALSTEIN® Refrigerators and Freezer**  
**Medical Refrigerator**  
**Model YR05083**



<b>Alarm</b>	High/low temperature	Y
	Power failure	Y
	Thermostat failure	Y
	Backup battery failure alarm	/
	Condenser fault alarm	/
	Ambient Temperatures	/
	Voltage abnormality alarm	/
	Door ajar	Y
	System abnormal alarm	/
	Power failure backup system(alarm)	8h
<b>Accessories</b>	Foot	Y
	Castor	/
	Outer door/Type	1/Glass/Low-E
	Door Lock	1
	Test hole QTY/diameter	1/25mm
	Lamp	LED
	USB port	/
	Temperature recorder	Optional
	Remote alarm port	Y
	RS232/485 port	/
<b>Other</b>	Certification	ISO9001/13485/14
		001;CE



#Letsgivemore

Thanks to your purchase, a donation will be made to a non-profit foundation that fights against breast cancer and helps most vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,  
 2 Rue Jean Lantier • 75001 Paris •  
 +33 1 78 95 87 89 / +33 6 80 76 07 10 •  
<https://kalstein.eu>  
**KALSTEIN FRANCE, S. A. S**





 **Kalstein**<sup>®</sup>  
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE  
ACQUISITION  
OF A KALSTEIN  
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación  
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



**Humatem**

Solidarité et coopération  
biomédicale


<https://www.humatem.org/>



**MANIAPURE  
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

**Articles Kalstein**  
(Social Networks and Official Websites)

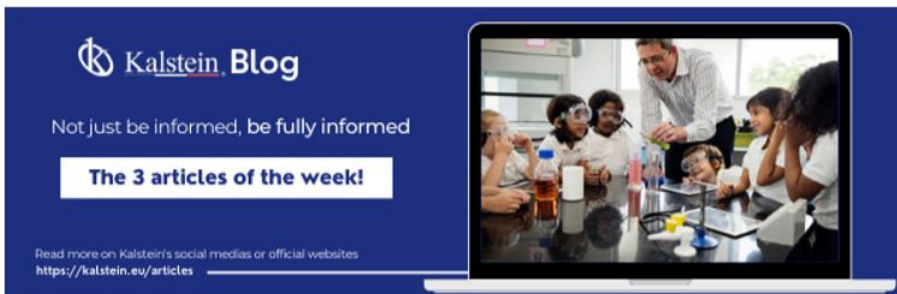

**Kalstein Blog**  
Not just be informed, be fully informed

**The 3 articles of the week!**

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**1**  
Read this ARTICLE NOW

**How does the climate affect your microscope?**

**Kalstein Blog**

**Article of the week:**  
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>



**2**  
Read this ARTICLE NOW

**What is the importance of the laboratory refrigerator?**

**Kalstein Blog**

**Article of the week:**  
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites  
<https://kalstein.eu/articles>

**Kalstein Training Videos**  
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Allianc...

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago



PARIS - FRANCE  
**HEADQUARTER**

2 Rue Jean Lantier  
75001 Paris - France

**Tlf:** +33 (0) 01 78 95 87 89

**Cel:** +33 (0) 6 80 76 07 10

[sales@kalstein.eu](mailto:sales@kalstein.eu)

<https://kalstein.eu/>  
<https://www.kalstein.fr/>

