




LABORATORY LINES SHAKING WATER BATH



Kalstein®

A different accompaniment, at your service

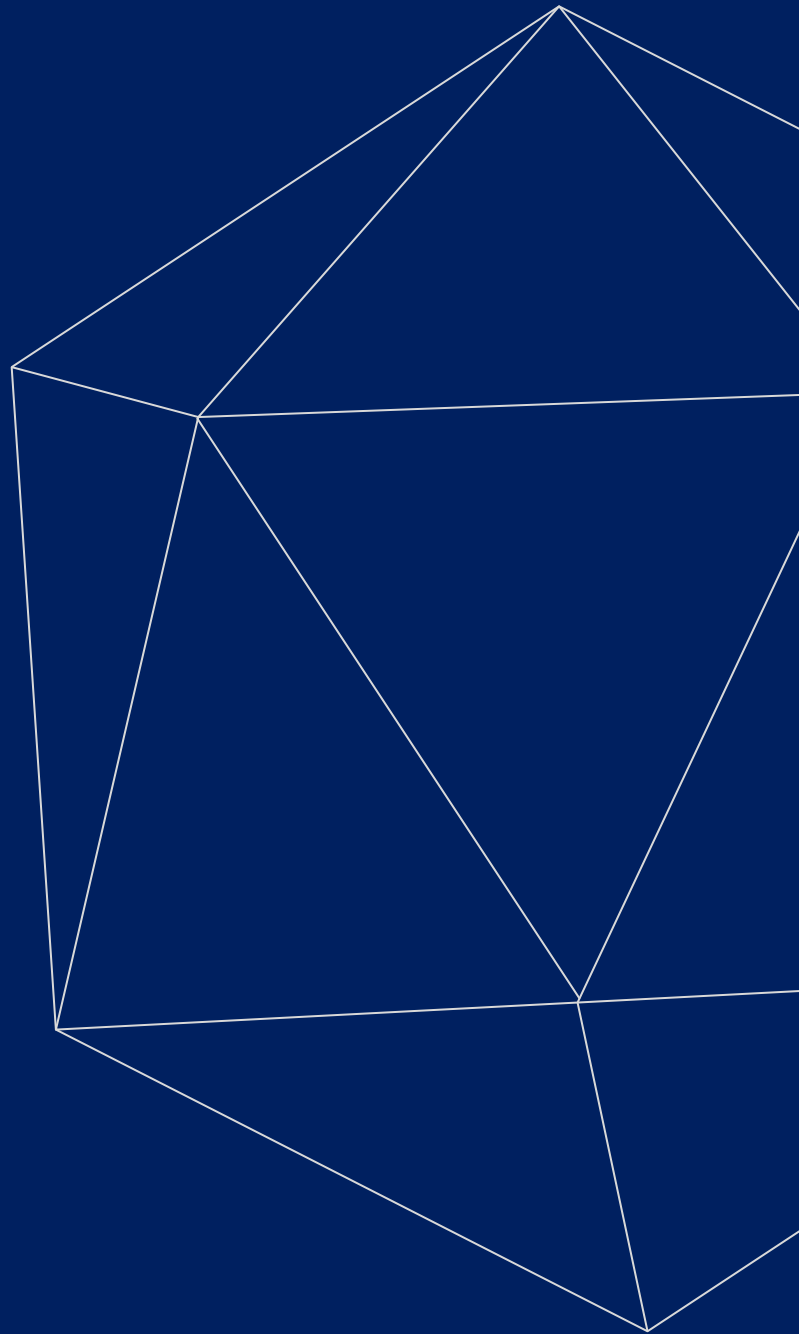



Kalstein®

A different accompaniment, at your service

OUR SERVICES

Benefits and Support

In Kalstein France, we take care of the full satisfaction of our customers, that is why we provide value-added services of the highest level based on our experience.



Online Inductions and Trainings

In any part of the world, receive your induction or training from our specialized team of engineers



Quick Response

Our work team is always available to response all your consults or questions, in order to support you in any situation.



#Letsgivemore ♥

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



Technical Support

Enjoy of personalized advice for the correct preventive and corrective maintenance of your equipment, thanks to Kalstein's manuals and articles, special catalogues and video tutorials.



Delivery Logistics

We take care of all the necessary logistics for the dispatch of your goods, whether is by sea, land or air.



Kalstein Worldwide

With more than 25 years growing with our customers, Kalstein's multiformat and modern content, is now present in more than 10 countries and increasing.



KALSTEIN® Water Baths

Shaking Water Bath

Model YR030L

Kalstein's shaking water bath applies to various large and medium-sized laboratories, medical, petrochemical industry, scientific research departments, biochemical, among others. It applies to various kinds of liquids as well as solid compounds and shaking culture.



Features:

- ✓ Ultra-silence design keeps it working quiet.
- ✓ Brushless motor cuts down the noise, lengthen the service life, and stabilize the torque output without maintenance.
- ✓ Integrated LCD controller, high accuracy sensor and precise temperature control
- ✓ Over temperature alarm prevents experiment temperature from exceeding the settings.
- ✓ Convenient PID Control makes it run with data set and time set, automatic changes of rotation direction, automatic shutdown, automatic date backups for power outage, quick start and stop, and motor shaft lock protection.



#Let'sgive more

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S





KALSTEIN® Water Baths Shaking Water Bath Model YR030L

Function	Operating Temp Range	RT~100°C
	Tem. Resolution Ratio	0.1°C
	Tem. Motion	±0.1°C
	Tem. Uniformity	±0.2°C
	Reciprocal Oscillation Amplitude	16 or 24mm (Ex-factory is 16mm)
	Reciprocating Frequency Range	20~180 (r/min)
	Reciprocating Frequency Accuracy	±1rpm
Structure	Water Tank	Mirror stainless steel
	Outer Shell	Cold rolling steel electrostatic spraying exterior
	Upper Cover	Mirror stainless steel
	Insulating Layer	Polyurethane
	Heater	Stainless Steel Heater
	Driving Mode	Crankshaft+ double link + four wheels
	Shelf	Mirror stainless steel
	Power rating	1.5kW
	Controller	Tem. Control Mode
Reciprocating Control Mode		
Setting Mode		Touch button setting
Tem. Display Mode		LCD Display
Speed Display Mode		
Timer		0~9999 min
Operation Function		Fixed value operation, timing function, auto stop, quickly stop start
Temperature Sensor		Pt100
Speed Sensor		Hoare
Additional Function	Deviation correction, menu keys locked, power-off memory	
Safety device	over temperature alarm, menu locked, water full or shortage protection, Motor shaft lock protection	
Specification	Water Tank Size (W*L*H)(mm)	500*310*200
	Exterior Size (W*L*H)(mm)	828*360*425
	Packing Size (W*L*H)(mm)	908*440*505
	Spring Shock Basket Size (mm)	395*250
	Shock Basket Largest Single Size (ml /piece)	100ml*12/250*8/500*6
	Inside Volume	31L
	Shock Basket Bearing	5kg/layer
	Current Rating (50Hz)	AC220V/6.8A
NW/GW (kg)	20/30	
Accessories	Water bath flask seat, tube rack	



#Letsgivemore

Thanks to your purchase, a donation is made to a non-profit foundation that fights against cancer and helps vulnerable communities.



All rights reserved © KALSTEIN France S. A. S.,
Optimum Business Center 450 Rue Baden Powell,
• 34000 Montpellier, France.
Tél: +33 467158849 / +33 680760710/ +33 663810023
<https://kalstein.eu>
KALSTEIN FRANCE, S.A. S





 **Kalstein**[®]
A different accompaniment, at your service

#LETSGIVEMORE 

**WITH THE
ACQUISITION
OF A KALSTEIN
EQUIPMENT**

YOU MAKE A CONTRIBUTION TO:



**Fundación
Jacinto Convit**

<https://www.jacintoconvit.org/>



<https://onetreepanted.org>



Humatem

Solidarité et coopération
biomédicale


<https://www.humatem.org/>



**MANIAPURE
FOUNDATION USA**

<https://www.maniapurefoundation.org/>

Articles Kalstein
(Social Networks and Official Websites)

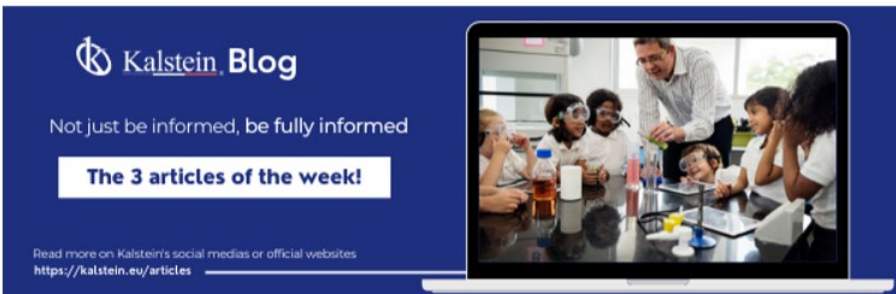

Kalstein Blog
Not just be informed, be fully informed

The 3 articles of the week!

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



PARIS - FRANCE
HEADQUARTER
2 Rue Jean Lantier
75001 Paris - France

Tlf: +33 (0) 01 78 95 87 89
Cel: +33 (0) 6 80 76 07 10
sales@kalstein.eu

<https://kalstein.eu/>
<https://www.kalstein.fr/>



1
Read this ARTICLE NOW

How does the climate affect your microscope?

microscope is an instrument that allows observing and measuring objects that are too small to be seen with the naked eye. This is achieved through a complex system of lenses and mirrors, which form

Kalstein Blog

Article of the week:
How does the climate affect your microscope?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>



2
Read this ARTICLE NOW

What is the importance of the laboratory refrigerator?

refrigerator is a device that is one of the most common pieces of laboratory equipment. Its purpose is to maintain, in a controlled environment (refrigerated space), various fluids and substances

Kalstein Blog

Article of the week:
What is the importance of the laboratory refrigerator?

Read more on Kalstein's social medias or official websites
<https://kalstein.eu/articles>

Kalstein Training Videos
(Youtube Channel and Official Websites)



Kalstein France - Jacinto Convit Foundation / Alliances

23 views • 3 years ago



Kalstein's Microscopes

445 views • 1 week ago

