

-86 °C Ultra Low Temperature Freezer



AUTO-CASCADE SYSTEM

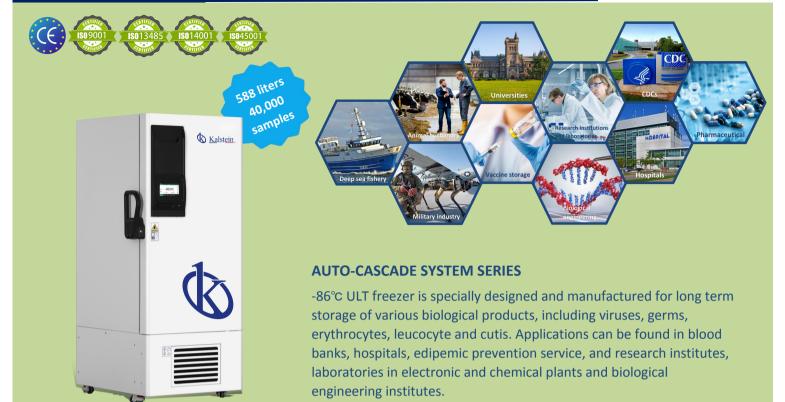

Pull down test at ambient temp.

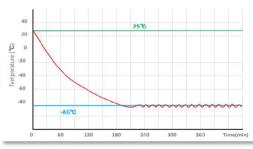

YR05110


Advanced thermal insulation

The insulation are designed with the latest generation of VIP technology to achieve thermal conductivity of 0.005 W(m.K). That makes ULTs more efficient in Insulation and less heat loss, so as to achieve the purpose of energy saving as well as space-saving.



Super-cooling technology


The powerful cooling system drives the unique mixed refrigeration system for fast cooling down. This system is stable, energy saving and easy maintenance.

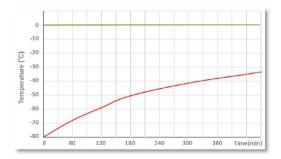

Pull up test at ambient temp. 25°C



EXTERNAL STRUCTURE

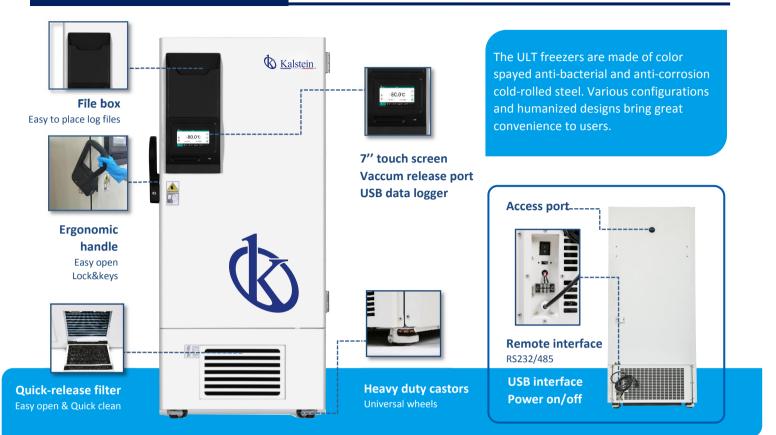


FEATURERS



7" LCD touch screen Easy operation and observation

Smart Operating System Authorities, alarms, curve & events

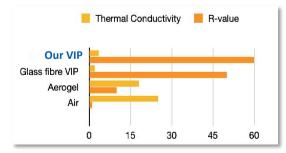





Energy saving More efficient cooling system





Advanced insulation tech Fumed Silicon V.I.P & PURF



Multi-layer sealing design Foamed inner doors with silicon gaskets



INTERNAL STRUCTURE



SPECIFICATIONS

Basic		Refrigeration System		Alarms	
Cooling method	Direct Cooling	Refrigerant	R4001(mixed gas)	High/low temperature	Y
Defrost	Manual	Temp. controller	Microprocessor	Power failure	Y
Temperature range	-40°C~-86°C	System	Auto-cascade	Abnormal voltage	Υ
Ambient temperature	10~32 ℃	Sensor	PTC	Sensor error	Υ
Display	7" LCD touch screen	Compressor brand/QTY	SECOP/2	Door ajar	Y
Access port QTY/Dia.	1/25mm	Capacity(L/cu.ft)	728/25.7	Filter screen check	Y
USB Intereface	Yes/1	Insulation	PURF+VIP	Thermostat failure	Υ
Remote alarm Access	RS232/485	Inner door material/QTY	Foamed SUS304/2	Low battery	Y
Outside material	Color sprayed steel	Racks	18(5x4)+6(6x4)	Condenser clean	Y
Inside material	SUS304	Boxes (2 Inches)	504	High ambient temperature	Y
Castor	Universal	Samples(2ml)	50,400	Power failure backup (alarm)	72h

Electrical Data		Dimensions		Accessories	
Voltage frequency(V/Hz)	220V 50Hz/60Hz	Package size(W*D*H)(mm)	1320*970*2185	Chart recorder	Optional
Power(W)	1400W	Exterior size(W*D*H)(mm)	1229*1012*1994.5	Data logger	Optional
Current(A)	10.37A	Interior size(W*D*H)(mm)	963*571*1350	CO2 back up system	Optional
Power consumption	23.44KWh/24h	Shipping CBM	2.79	LN2 back up system	Optional
Noise(db)	53	N.W/G.W(kg)	374/394	Gloves	Optional

ACCESSORIES









LN2/CO2 backup system

Data logger

Chart recorder

Racks & Boxes

Gloves